


# ÜBERSICHT BOHRPROTOKOLL ELEMENT

PF 3.0 enossal Ø 2.5	PF Ø 3.5 enossal Ø 3.5	PF Ø 4.0 enossal Ø 4.0	PF Ø 4.5 enossal Ø 4.2	PF Ø 5.0 enossal Ø 5.0	PF Ø 6.0 enossal Ø 6.0	Drehzahl U/min
2.0 Pilotbohrer	2.0 Pilotbohrer	2.0 Pilotbohrer	2.0 Pilotbohrer	2.0 Pilotbohrer	2.0 Pilotbohrer	800
 Längenspezifische Profilbohrer	2.8 VECTOdrill	2.8 VECTOdrill	2.8 VECTOdrill	2.8 VECTOdrill	2.8 VECTOdrill	600
		3.5 VECTOdrill	3.5 VECTOdrill	3.5 VECTOdrill	3.5 VECTOdrill	500
					4.3 VECTOdrill	4.3 VECTOdrill
5.3 VECTOdrill	5.3 VECTOdrill				300	
	3.5/3.5 ELEMENT Profilbohrer*	4.0/4.0 ELEMENT Profilbohrer*	4.2/4.5 ELEMENT Profilbohrer*	5.0/5.0 ELEMENT Profilbohrer*	6.0/6.0 ELEMENT Profilbohrer*	250-300

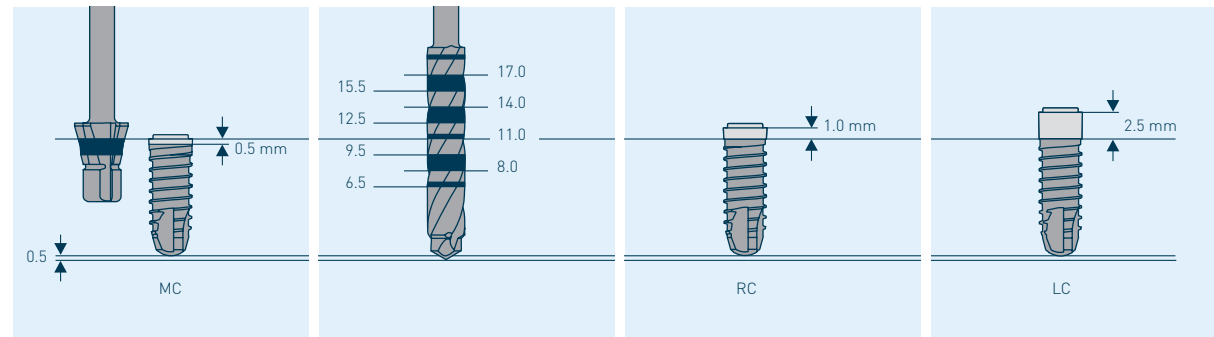

\* nur bei ELEMENT MC